

Product Information

N6,2'-O-Dibutyryl adenosine 3',5'-cyclic monophosphate sodium salt

≥96% (HPLC), powder

D0627

Product Description


CAS Registry Number: 16980-89-5

Molecular Formula: C₁₈H₂₃N₅O₈PNa

Formula Weight: 491.37 (anhydrous)

Synonyms: Dibutyryl cAMP sodium salt, Bucladesine sodium salt, Dibutyryl cyclic-AMP sodium salt, Bucladesine, Dibutyryl cAMP

Structure:

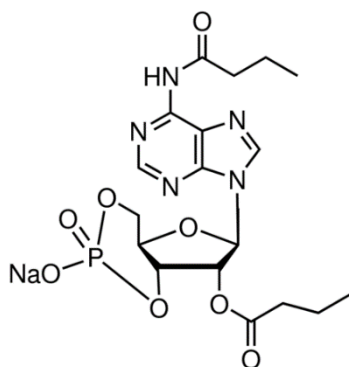

Dibutyryl cAMP is an analog of cAMP (cyclic AMP; adenosine 3',5'-cyclic monophosphate) that mimics the action of endogenous cAMP.^{1,2} Compared to cAMP, the lipophilic nature of dibutyryl cAMP gives it greater cell permeability, and greater resistance to hydrolysis by cAMP phosphodiesterases.^{3,4} Known to activate cAMP-dependent protein kinases and to inhibit phosphodiesterases, dibutyryl cAMP is used to probe signal transduction pathways.⁵ This product is synthetically prepared and chromatographically purified.¹

Dibutyryl cAMP is widely used in cell culture,⁶⁻⁹ such as for mediation of cell differentiation. Several publications^{6,7,10-20} and dissertations^{8,9,21-26} have cited use of product D0627 in their research protocols.

Absorbance: 273 nm in 0.1 M phosphate buffer (pH 7.0)

E_{mM}^{273} : 16.6 (0.1 M Phosphate, pH 7.0)

A_{250}/A_{260} : 0.75

A_{280}/A_{260} : 1.15

Preparation Instructions

This product is soluble in water at 100 mg/mL. With reconstituted solutions, because the 2'-O-butyryl group hydrolyzes at pH ≥ 8.5, pH ≥ 8.5 solutions should be avoided.²⁷ While we have not tested solution stability on this reagent, several publications have indicated storage of stock solutions of dibutyryl cAMP at -20 °C.^{28,29}

Storage/Stability

Dibutyryl cAMP, as supplied, is sensitive to light and to moisture. It is recommended to store this product at -20 °C.

Precautions and Disclaimer

For R&D use only. Not for drug, household, or other uses. Please consult the Safety Data Sheet for information regarding hazards and safe handling practices.

References

1. Posternak, T., and Weimann, G., *Meth. Enzymol.*, **38**, 399-409 (1974).
2. *The Merck Index*, 11th ed., Entry #1448, 221 (1989).
3. Henion, W.F. *et al.*, *Biochem. Biophys. Acta*, **148(1)**, 106-113 (1967).
4. Swislocki, N.I., *Anal. Biochem.*, **38(1)**, 260-269 (1970).
5. Schwede, F. *et al.*, *Pharmacol. Ther.*, **87(2-3)**, 199-226 (2000).
6. Goorha, S., and Reiter, L.T., *Curr. Protoc. Hum. Genet.*, **92**, 21.6.1 – 21.6.10 (2017).

7. Zabolocki, M. *et al.*, *Nat. Commun.*, **11(1)**, 5550 (2020).
8. McGinnis, Lauren Ashley, "Proteins that mediate coordination of the cell cortex and the meiotic spindle in mammalian oocytes: insights into cellular factors affecting successful reproduction". Johns Hopkins University, Ph.D. dissertation, pp. 136, 199 (2016).
9. Kim, Jihyun, "Molecular mechanisms of myelin dysfunctions in the nervous system". Rutgers, The State University of New Jersey, Ph.D. dissertation, p. 33 (2017).
10. Klingler, W. *et al.*, *Anesthesiology*, **97(5)**, 1059-1066 (2002).
11. Akama, K.T., and McEwen, B.S., *J. Neurosci.*, **23(6)**, 2333-2339 (2003).
12. Huang, H.W. *et al.*, *Am. J. Respir. Cell Mol. Biol.*, **38(4)**, 473-482 (2008).
13. Abbassi, L. *et al.*, *Biol. Reprod.*, **94(5)**, 102 (2016).
14. Abe, T. *et al.*, *Int. J. Oncol.*, **54(2)**, 558-558 (2019).
15. Fitzgerald, H.C. *et al.*, *Proc. Nat. Acad. Sci. USA*, **116(46)**, 23132-23142 (2019).
16. Semenov, A.N. *et al.*, *Front. Physiol.*, **10**, 923 (2019).
17. Vajhøj, C. *et al.*, *PLoS One*, **16(12)**, e0261536 (2021).
18. Yates, P.L. *et al.*, *Cell. Mol. Life Sci.*, **78(21-22)**, 6941-6961 (2021).
19. Yin, X. *et al.*, *STAR Protoc.*, **2(2)**, 100405 (2021).
20. Raabe, F.J. *et al.*, *Cells*, **11(2)**, 241 (2022).
21. Polikov, Vadim Steven, "An *In Vitro* Model of the Brain Tissue Reaction to Chronically Implanted Recording Electrodes Reveals Essential Roles for Serum and bFGF in Glial Scarring". Duke University, Ph.D. dissertation, p. 98 (2009).
22. De, Santanu, "Protein 14-3-3 (YWHA) isoforms and their roles in regulating mouse oocyte maturation". Kent State University, Ph.D. dissertation, p. 134 (2014).
23. Mukherjee-Clavin, Bipasha, "Modeling Charcot Marie Tooth 1a with human pluripotent stem cells". Johns Hopkins University, Ph.D. dissertation, p. 134 (2015).
24. Abernathy, Daniel Gene, "Brain Enriched microRNAs Open the Neurogenic Potential of Adult Human Fibroblasts". Washington University, Ph.D. dissertation, p. 102 (2017).
25. Breitmeyer, Ricarda, "Neuron-microglia interactions induce aberrant inflammatory mechanisms in schizophrenia". Eberhard Karls Universität Tübingen, Dr. rer. nat. dissertation, p. 18 (2020).
26. R'Bibo, Léa Emma Sophie, "Investigating the Role of GTP Cyclohydrolase I Mutations and the Tetrahydrobiopterin Pathway in Parkinson's Disease". University College London, Ph.D. dissertation, p. 75 (2020).
27. Dawson, R.M.C. *et al.*, *Data for Biochemical Research*, 3rd edition. Clarendon Press (Oxford, UK), pp. 78-79 (1986).
28. Wilding, T.J. *et al.*, *J. Neurosci.*, **15(5)**, 4124-4132 (1995).
29. Major, T. *et al.*, *Curr. Protoc. Stem Cell Biol.*, **39**, 1H.10.1 – 1H.10.23 (2017).

Notice

We provide information and advice to our customers on application technologies and regulatory matters to the best of our knowledge and ability, but without obligation or liability. Existing laws and regulations are to be observed in all cases by our customers. This also applies in respect to any rights of third parties. Our information and advice do not relieve our customers of their own responsibility for checking the suitability of our products for the envisaged purpose.

The information in this document is subject to change without notice and should not be construed as a commitment by the manufacturing or selling entity, or an affiliate. We assume no responsibility for any errors that may appear in this document.

Technical Assistance

Visit the tech service page at SigmaAldrich.com/techservice.

Standard Warranty

The applicable warranty for the products listed in this publication may be found at SigmaAldrich.com/terms.

Contact Information

For the location of the office nearest you, go to SigmaAldrich.com/offices.

The life science business of Merck KGaA, Darmstadt, Germany operates as MilliporeSigma in the U.S. and Canada.

MilliporeSigma, and Sigma-Aldrich are trademarks of Merck KGaA, Darmstadt, Germany or its affiliates. All other trademarks are the property of their respective owners. Detailed information on trademarks is available via publicly accessible resources.

© 2022 Merck KGaA, Darmstadt, Germany and/or its affiliates. All Rights Reserved.

**MILLIPORE
SIGMA**