

RABBIT ANTI-CYP450 Pan 2C  
POLYCLONAL ANTIBODY


CATALOG NUMBER: AB10324 QUANTITY: 100 µL

LOT NUMBER: HOST/ISOTYPE: Rabbit

IMMUNOGEN: Recombinant CYP450 2C19

APPLICATIONS: Western Blot Analysis: Western Blot Analysis of AB10324 at a dilution of 1:5,000.

Lane 1 Recombinant CYP450 2C8  
Lane 2 Recombinant CYP450 2C9  
Lane 3 Recombinant CYP450 2C18  
Lane 4 Recombinant CYP450 2C19  
Arrow indicates expected band at 56 kDa

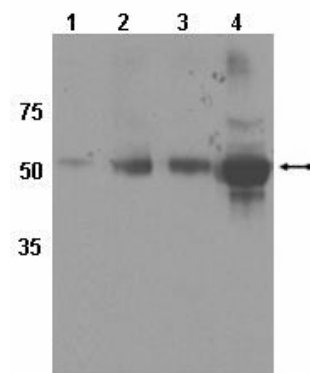

Optimal working dilutions must be determined by the end user.

SPECIES REACTIVITY: Human

FORMAT: Raw Rabbit Serum

PRESENTATION: Raw Rabbit Serum with 0.1% Sodium Azide

STORAGE/HANDLING: Maintain at 2-8°C for up to 12 months from date of receipt.

REFERENCES: Peng, Fu-Chuo, *et al* (2006). Territrems B and C metabolism in human liver microsomes: major role of CYP3A4 and CYP3A5. *Toxicology* **218**: 172-85Ko, J W, *et al* (2000). In vitro inhibition of the cytochrome P450 (CYP450) system by the antiplatelet drug ticlopidine: potent effect on CYP2C19 and CYP2D6. *British journal of clinical pharmacology* **49**: 343-51.**Important Note:**

*During shipment, small volumes of product will occasionally become entrapped in the seal of the product vial. For products with volumes of 200 µL or less, we recommend gently tapping the vial on a hard surface or briefly centrifuging the vial in a tabletop centrifuge to dislodge any liquid in the container's cap.*

**For research use only; not for use as a diagnostic.**