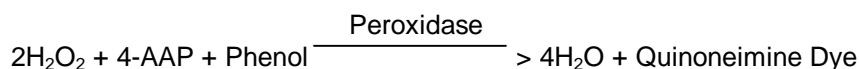
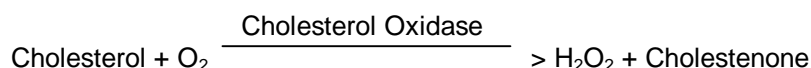


SIGMA QUALITY CONTROL TEST PROCEDURE
Enzymatic Assay of CHOLESTEROL ESTERASE
(EC 3.1.1.13)

PRINCIPLE:

Abbreviation:

4-AAP = 4-Aminoantipyrine

CONDITIONS: T = 37°C, pH = 7.0, A_{500nm}, Light path = 1 cm

METHOD: Continuous Spectrophotometric Rate Determination

REAGENTS:

- A. 400 mM Potassium Phosphate Buffer, pH 7.0 at 37°C
(Prepare 100 ml in deionized water using Potassium Phosphate, Monobasic, Anhydrous, Sigma Prod. No. P-5379. Adjust to pH 7.0 at 37°C with 1 M KOH.)
- B. 0.9% (w/v) Sodium Chloride Solution
(Prepare 25 ml in deionized water using Sodium Chloride, Sigma Prod. No. S-9625.)
- C. 8.6 mM Cholesteryl Oleate Solution (Chol-Oleate)
(Prepare 10 ml by first dissolving the Cholesteryl Oleate, Sigma Prod. No. C-9253, in 1 ml of Polyoxyethylene, 9 Lauryl Ether, Sigma Prod. No. P-9641. Stir gently, with heat, until the solution is clear and colorless. Then add 9 ml of hot Reagent B and continue for 5 minutes. Allow the solution to return to ambient temperature prior to use. The solution clears upon cooling.)

**Enzymatic Assay of CHOLESTEROL ESTERASE
(EC 3.1.1.13)**

REAGENTS: (continued)

- D. 15% (w/v) Taurocholic Acid Solution (Tauro)
(Prepare 10 ml in deionized water using Taurocholic Acid, Sodium Salt, Sigma Prod. No. T-4009.)
- E. 15% (w/v) Cholic Acid Solution (Chol)
(Prepare 10 ml in deionized water using Cholic Acid, Sodium Salt, Sigma Prod. No. C-1254.)
- F. 1.76% (w/v) 4-Aminoantipyrine Solution (4-AAP)
(Prepare 1 ml in deionized water using 4-Aminoantipyrine, Free Base, Sigma Prod. No. A-4382.)
- G. 5% (w/v) Phenol Solution (Phenol)
(Prepare 10 ml in deionized water using Phenol, Sigma Prod. No. P-4161.)
- H. Cholesterol Oxidase Enzyme Solution (Chol Oxid)
(Immediately before use, prepare a solution containing 20 - 30 units/ml of Cholesterol Oxidase, Sigma Prod. No. C-1512, in cold Reagent A.)
- I. Peroxidase Enzyme Solution (POD)
(Immediately before use, prepare a solution containing 40 - 60 units/ml of Peroxidase Type II from Horseradish, Sigma Prod. No. P-8250 in cold deionized water.)
- J. Cholesterol Esterase Enzyme Solution
(Immediately before use, prepare a solution containing 0.25 - 0.85 unit/ml of Cholesterol Esterase in cold Reagent A.)

PROCEDURE:

Pipette (in milliliters) the following reagents into suitable cuvettes¹:

	<u>Test</u>	<u>Blank</u>
Reagent A (Buffer)	2.10	2.10
Reagent D (Tauro)	0.05	0.05
Reagent E (Chol)	0.05	0.05
Reagent I (POD)	0.10	0.10
Reagent C (Chol-Oleate)	0.50	0.50
Reagent G (Phenol)	0.05	0.05

**Enzymatic Assay of CHOLESTEROL ESTERASE
(EC 3.1.1.13)**

PROCEDURES: (continued)

Mix by inversion and equilibrate to 37°C. Then add:

Reagent F (4-AAP)	0.05	0.05
Reagent H (Chol Oxid)	0.05	0.05

Mix by inversion and obtain the baseline at 500 nm. After approximately 5 minutes add:

	<u>Test</u>	<u>Blank</u>
Reagent J (Cholesterol Esterase)	0.05	-----
Reagent A (Buffer)	-----	0.05

Immediately mix by inversion and record the increase in A_{500nm} for approximately 10 minutes.²
Obtain the $\Delta A_{500nm}/\text{minute}$ using the maximum linear rate for both the Test and Blank.

CALCULATIONS:

$$\text{Units/ml enzyme} = \frac{(\Delta A_{500nm}/\text{min Test} - \Delta A_{500nm}/\text{min Blank})(3)(df)}{(0.5)(13.78)(0.05)} >$$

3 = Total volume (in milliliters) of assay

df = Dilution factor

0.5 = Conversion factor based on one mole of H_2O_2 produces half a mole of Quinoneimine Dye

0.05 = Volume (in milliliters) of enzyme used

13.78 = Millimolar extinction coefficient of Quinoneimine Dye at 500 nm under the assay conditions

$$\text{Units/mg solid} = \frac{\text{units/ml enzyme}}{\text{mg solid/ml enzyme}}$$

$$\text{Units/mg protein} = \frac{\text{units/ml enzyme}}{\text{mg protein/ml enzyme}}$$

UNIT DEFINITION:

One unit will hydrolyze 1.0 μmole of cholesteryl oleate to cholesterol and oleic acid per minute at pH 7.0 at 37°C in the presence of taurocholate.

**Enzymatic Assay of CHOLESTEROL ESTERASE
(EC 3.1.1.13)**

FINAL ASSAY CONCENTRATION:

In a 3.00 ml reaction mix, the final concentrations are 287 mM potassium phosphate, 0.25% (w/v) taurocholic acid, 0.25% (w/v) cholic acid, 4 - 6 units peroxidase, 1.4 mM cholesteryl oleate, 1.7% (v/v) polyoxyethylene 9 lauryl ether, 0.14% (w/v) sodium chloride, 0.083% (w/v) phenol, 0.03% (w/v) 4-aminoantipyrine, 1 - 1.5 units cholesterol oxidase and 0.013 - 0.043 unit cholesterol esterase.

REFERENCE:

Allain, C.C. et al., (1974) *Clinical Chemistry*, **20**, 470-475

NOTES:

1. Add the reagents in the order written.
2. The fastest rate is usually between 4-8 minutes after addition of the Cholesterol Esterase.
3. Cholesterol Oxidase Unit Definition: One unit will convert 1.0 μ mole of cholesterol to 4-cholesten-3-one per minute at pH 7.5 at 25°C.
4. Peroxidase Unit Definition: One unit will form 1.0 mg purpurogallin from pyrogallol in 20 sec at pH 6.0 at 20°C.
5. This assay is based on the cited reference.
6. Where Sigma Product or Stock numbers are specified, equivalent reagents may be substituted.

Sigma warrants that the above procedure information is currently utilized at Sigma and that all Sigma-Aldrich, Inc. products conform to the information in this and other Sigma-Aldrich, Inc. publications. Purchaser must determine the suitability of the information and product(s) for their particular use. Additional terms and conditions may apply. Please see reverse side of the invoice or packaging slip.