

Product Information

3-Amino-9-ethylcarbazole

≥95% (HPLC), powder

A5754

Product Description


CAS Number: 132-32-1

Synonyms: AEC, 9-ethylcarbazol-3-amine

Molecular Weight: 210.27

Molecular Formula: C₁₄H₁₄N₂

Melting Point: 98.0 – 100.0 °C

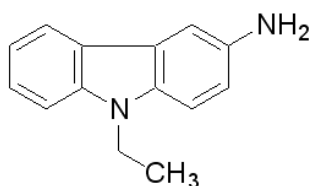

3-Amino-9-ethylcarbazole (AEC) is a chromogen and enzyme substrate that is suitable for use in immunoblotting and immunohistochemical staining procedures which utilize horseradish peroxidase (HRP) conjugates.¹⁻³ The reaction of AEC with HRP produces an insoluble end product that is red in color and can be observed visually, in 0.05 M acetate buffer at pH 5. An aqueous mounting medium should be used with this product, as the end product is alcohol-soluble.

AEC has been shown to work better than benzidine in myeloperoxidase staining¹ and to be less sensitive to oxidases than phenazine methosulfate and nitroblue tetrazolium.⁴ However, AEC is not as sensitive for HRP as 2,2'-azino-bis(3-ethylbenzthiazoline-6-sulfonic acid (ABTS) or *o*-phenylenediamine (OPD).⁵ AEC can also be used to detect cytochrome oxidase.⁶ AEC can be used in combination with *N*-phenyl-*p*-phenylenediamine using cobalt ion chelation to enhance the observed color.⁷

Several dissertations⁸⁻¹⁷ have cited use of product A5754 in their protocols.

Precautions and Disclaimer

For R&D use only. Not for drug, household, or other uses. Please consult the Safety Data Sheet for information regarding hazards and safe handling practices.

Storage/Stability

Store this product at 2-8 °C.

Solubility

This product is tested for solubility at 50 mg/mL in a 2:1 ethanol:acetic acid mixture. AEC is also soluble in DMF (*N,N*-dimethylformamide), at concentrations like 50 mg/mL.¹⁷ Another publication reports preparation of AEC stock solutions in DMSO at 20 mg/mL.¹⁸

Solutions of AEC in DMF may be stored at 4 °C when stored in the dark and kept tightly sealed. One publication reports storage of at 20 mg/mL AEC stock solutions in DMSO at -20 °C, protected from light by foil covering.¹⁸

Usage

For use as a peroxidase (HRP) substrate

- Dissolve 25 mg of AEC in 2.5 mL of DMF.
- This AEC solution is then added to 47.5 mL of 50 mM acetate buffer, pH 5.0 (final AEC concentration of 0.5 mg/mL).
- Immediately before use, add 25 µL of 30% hydrogen peroxide to this solution.

References

1. Kaplow, L.S., *Am. J. Clin. Path.*, **63(3)**, 451 (1975).
2. Burstone, M.S., and Weisburger, E.K., *J. Histochem. Cytochem.*, **9(4)**, 349-355 (1961).

3. Graham, R.C., Jr. *et al.*, *J. Histochem. Cytochem.*, **13(2)**, 150-152 (1965).
4. Feinstein, R.N., and Lindahl, R., *Anal. Biochem.*, **56(2)**, 353-360 (1973).
5. Al-Kaissi, E., and Mostratos, A., *J. Immunol. Meth.*, **58(1-2)**, 127-132 (1983).
6. Burstone, M.S., *Enzyme histochemistry and its application in the study of neoplasms*. Academic Press (New York, NY), pp. 429-443 (1962).
7. Burstone, M.S., *J. Histochem. Cytochem.*, **8(1)**, 63-70 (1960).
8. Irwin, Judith A., "Male competition and outcrossing rate in a hermaphrodite plant". University of St. Andrews, Ph.D. dissertation, Appendix A (1990).
9. Hemery, Gabriel E., "*Juglans regia* L: genetic variation and provenance performance". University of Oxford, Ph.D. dissertation, p. 170 (2000).
10. Fernando, Neluka, "Investigations into the variation of strains of *Helicobacter pylori* and the Host response to infection from different geographical locations". University College London, Ph.D. dissertation, p. 110 (2001).
11. Constantinescu, Ana Alina, "Degradation of the endothelial glycocalyx by atherogenic factors: Microvascular functional implications". Universiteit van Amsterdam, Ph.D. dissertation, p. 98 (2002).
12. Michel, Brady, "Early Peripheral Immunological Events Dictate Chronic Wasting Disease". Colorado State University, Ph.D. dissertation, p. 96 (2013).
13. Riederer, Isabelle Marie Mathilde, "Interaction of Endothelial Cells and Natural Killer Cells". Technischen Universität München, Dr. med. dissertation, p. 29 (2014).
14. Oudijk, Lindsey, "Immunohistochemical and Molecular Studies on the Pathogenesis of Pheochromocytomas and Paragangliomas". Erasmus Universiteit Rotterdam, Ph.D. dissertation, p. 126 (2015).
15. Elfaitwri, Tehani Ali Omer, "VHL and Hypoxia regulated genes in sympatho-adrenal cell lineage Development and Disease". Albert-Ludwigs-Universität Freiburg im Breisgau, Dr. med. dissertation, p. 21 (2016).
16. Budd, Emma, "MicroRNAs in Osteoarthritis and Chondrogenesis". University of Southampton, Ph.D. dissertation, pp. 90, 92 (2016).
17. Câmara, Nakita, "Stress Cardiomyopathy in Cetaceans: Histological, Histochemical, Immunohistochemical and Biochemical Studies". Universidad de Las Palmas de Gran Canaria, Ph.D. dissertation, p. 121 (2020).
18. Nirenberg, S., *Methods Mol. Biol.*, **135**, 475-480 (2000).

Notice

We provide information and advice to our customers on application technologies and regulatory matters to the best of our knowledge and ability, but without obligation or liability. Existing laws and regulations are to be observed in all cases by our customers. This also applies in respect to any rights of third parties. Our information and advice do not relieve our customers of their own responsibility for checking the suitability of our products for the envisaged purpose.

The information in this document is subject to change without notice and should not be construed as a commitment by the manufacturing or selling entity, or an affiliate. We assume no responsibility for any errors that may appear in this document.

Technical Assistance

Visit the tech service page at [SigmaAldrich.com/techservice](https://www.sigmaaldrich.com/techservice).

Terms and Conditions of Sale

Warranty, use restrictions, and other conditions of sale may be found at [SigmaAldrich.com/terms](https://www.sigmaaldrich.com/terms).

Contact Information

For the location of the office nearest you, go to [SigmaAldrich.com/offices](https://www.sigmaaldrich.com/offices).

The life science business of Merck KGaA, Darmstadt, Germany operates as MilliporeSigma in the U.S. and Canada.

MilliporeSigma, and Sigma-Aldrich are trademarks of Merck KGaA, Darmstadt, Germany or its affiliates. All other trademarks are the property of their respective owners. Detailed information on trademarks is available via publicly accessible resources.
© 2022 Merck KGaA, Darmstadt, Germany and/or its affiliates. All Rights Reserved.

A5754pis Rev 07/22 CMH,RXR,GCY

**MILLIPORE
SIGMA**